


SinoMRI series

| MRI Injectors

SinoMRI

NEW

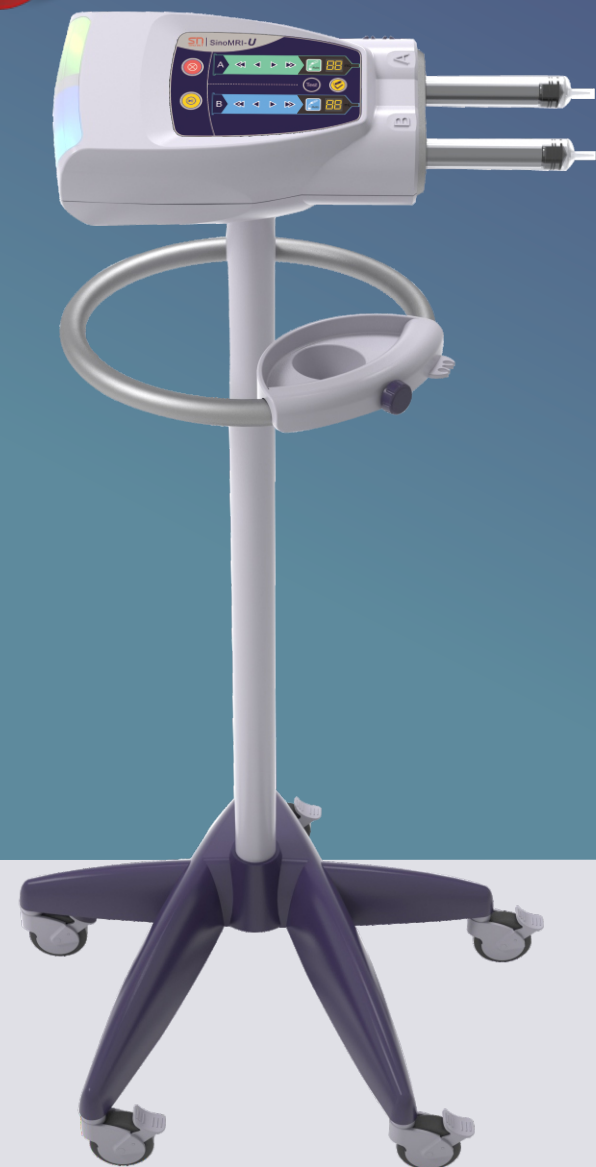


SinoMRI-U



Unique *display* function

- Injection head digital display the volume of contrast media

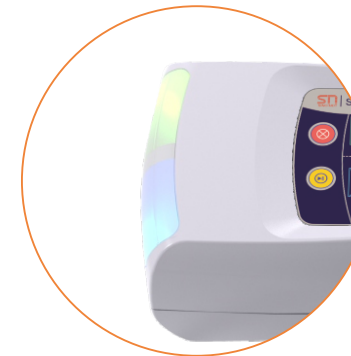

LED working status light

- LED shows dark green or dark blue, means the system is in standby state
- LED shows high light green or high light blue, means the system is in ready state
- LED shows blink high light green or blink high light blue, means the system is in injection state
- LED shows blink red, means the system is in alarm



Function

1. "Direct pressure sensor" real-time detection of injecting pressure, Automatic alarm and stop injection when abnormal pressure.
2. "Air exhaustion confirmation" to make sure injection after exhausting air.
3. 12.1Inch color touch screen, Graphics User Interface.
4. LED working status light and LED indicator light specially shows standby, ready, injection and alarm state.
5. It is convenient to use soft or hard "stop" buttons to stop injection.
6. The anti-leakage design of injector head can prevent injector head being damaged by liquid.